



ENTONOX CD CYLINDER USAGE

1

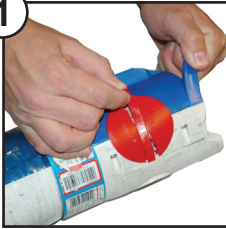

Remove plastic tabs.

2

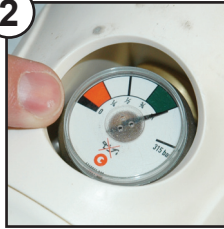

Check cylinder contents.

3




If temperature is below 5°C, activate the warming pack and place around the cylinder. Put cylinder and pack back inside the barrel bag and re-seal.

4



Regardless of temperature, invert cylinder a few times.

5

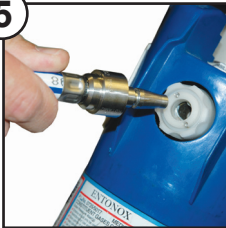


Insert bayonet hose.

6



Turn on main valve.

7

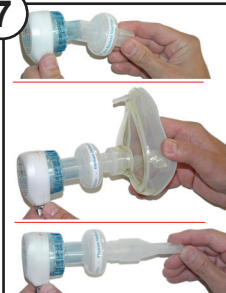

Attach the filter to the regulator.

Attach the mask OR mouthpiece to the filter.

8



If possible, patient self administers deep slow breaths.

NOTES:

- Pain relief occurs between 1 and 10 minutes and effect clears in approx. 2 minutes.
- Heat pack can be recharged by simmering in boiling water for 15 minutes.
- The filter is for single patient use and must be disposed of after treatment.
- Using the filter negates the requirement for cleaning the patient demand valve.
- Mouthpiece is single patient use and disposable, face mask requires cleaning after use.